


SKETCH

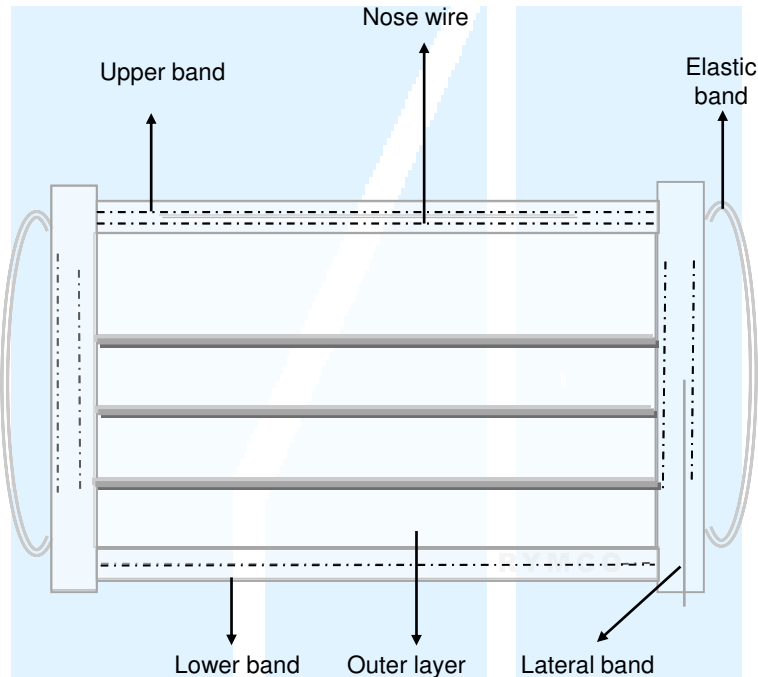

**STANDARD EARLOOP
FACE MASK BFE 98%**

Page
1/2

DEVICE DESCRIPTION:

Non-esterile disposable face mask for single use.

MAIN COMPONENTS

outer layer: Made of nonwoven fabric or cellulose with a minimum grammage of 13 gsm.

Middle layer - Barrier: : Made of nonwoven fabric filter with a minimum grammage of 10 gsm. (Bacterial filtration efficiency $\geq 98\%$)

Inner layer: Made of nonwoven fabric or cellulose with a minimum grammage of 17 gsm.

Lateral bands: Made of nonwoven fabric with a minimum grammage of 20 gsm.

Nose wire: Made of plastic coated wire (Polyethylene), 3.0 mm wide and wire gauge of 0.5mm.

Elastic band: Made with Spandex and Polyester, does not irritate the skin.

Sealed with ultrasound

Available colors: Yellow, Blue, White, Lilac, Green.

Face mask Dimensions:

Mask length: 172 mm -178 mm

Mask width: 98 mm -102 mm

Wire length: minimum 95 mm

Elastic length: 140 mm – 180 mm

lateral band width: 8 mm – 12 mm.

PACKING UNITS

25 units / Folding box

50 units / Folding box

20 Folding box / Corrugated box

1000 units / Corrugated box

Plastic bag x 32 units, Corrugated box x 50 bags, Corrugated box x 1600 units.

Double-sided printed bag x 1 unit, corrugated box x 1000 units.

Double-sided printed bag x 5 units, corrugated box x 1000 units.

Double-sided printed bag x 10 units, corrugated box x 1000 units.

Double-sided printed bag x 25 units, corrugated box x 1000 units.

Double-sided printed bag x 50 units, corrugated box x 1000 units.

VISUAL INSPECTION

Absence of defects: Facial masks must be free of stains, foreign particles and strands. Facial masks must be free of holes. According to internal manufacturing standards.

PHYSICAL TESTS

Dimensions: mask length, mask width, elastic length, wire length and lateral band width. According to internal manufacturing standards.

Mechanical resistance: determines the safety of the face mask joints when subjected to external tension forces. According to internal manufacturing standards.

INDICATIONS

- Use the Facial Mask according to the requirements of the health institution.
- Use aseptic technique and open the packaging.
- Discard once it has been used, according to the institution's protocols.


USES AND PRECAUTIONS

The Facial Mask for single use is a device used to reduce the risk of contact with polluting agents or risk to the health of medical personnel and patients, as well as to avoid contamination with fluids.

The Facial Mask is used by doctors, professional nurses, nursing assistants and other personnel that require it. If the Facial Mask is not in good condition, do not use this product. This element should only be used once, discard the facial mask once it has been used.

**STANDARD EARLOOP
FACE MASK BFE 98%**

Page
2/2

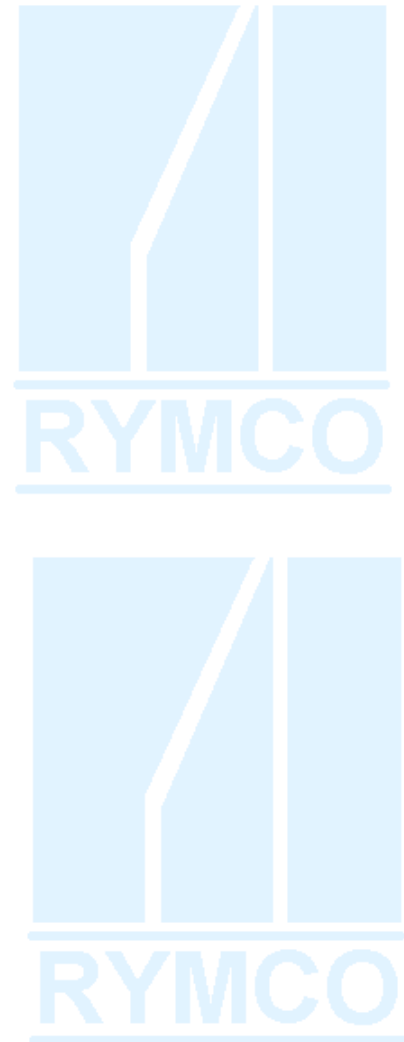

Once downloaded or printed this document will be considered an uncontrolled copy. Verify its version and validity before using. The original document is in the custody of the Documentation Administrator.
Document owned by RYMCO S.A. total or partial reproduction is prohibited.

DATE OF ISSUE
2020 - 02 - 10

Barranquilla - Colombia
Calle 80 No. 78B - 51 A.A. 3130
www.rymco.com.co
E-mail: calidad@rymco.com.co
Quality Department

PRODUCED BY:
Cristian H.
CRISTIAN HERNANDEZ D.
DEVELOPMENT AND SPECIFICATION
ANALYST

REVIEW:
Daniel Cudriz Z.
DANIEL CUDRIZ Z.
QUALITY CONTROL COORDINATOR

APPROVAL:
Carlos Solano G.
CARLOS SOLANO G.
QUALITY AND DEVELOPMENT DIRECTOR